Brugsanvisning

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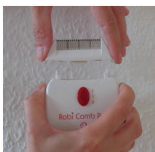
(English)

ROBI COMB & ROBI COMB PRO

- A product that detects, eliminates and repels headlice without the use of chemicals.

Congratulations on your decision to join the growing number of ROBI COMB PRO purchasers using this innovative method to control lice without chemicals. This new generation electronic lice-eliminator, identifies and kills lice when in contact with the combining unit teeth of the comb. Thanks to ROBI COMB PRO's patented design, only the lice are harmed and it is completely safe to use for both adults and children aged 3 years or older.

How does ROBI COMB PRO work? Liceguard ROBI COMB PRO runs on AA battery and consists of a comb with metal teeth. When headlice are detected between the two teeth, it gives off a small electric current, which passes through the lice and kills it. The teeth on the comb are protected, so that the electric current will not be felt on the head.

ROBI COMB PRO had been sold in millions, in over 30 countries throughout the world.

For further treatment of headlice we recommend LiceGuard's other products: LiceGuard Shampoo, Liceegg comb, Lice Repellent Spray.

Guarantee: LiceGuard makes every effort to ensure that it's products are manufactured without any defect. However in the rare event of defects due to faulty workmanship or materials, we guarantee to repair or replace the product, or any part thereof, free of charge.

This guarantee is granted to the original purchaser only and it is valid for 2 years from the date of purchase. This guarantee does not apply if the product is altered, repaired or otherwise tampered with by unauthorized persons and/ or is used other than in accordance with the manufacturers instructions and/ or is subjected to use other than reasonable use.

Should a defect occur please contact your place of purchase along with the original receipt of purchase. Follow the instruction of use.

Instructions for use of ROBI COMB & ROBI COMB PRO

1. Insert a AA 1.5 Volt battery (not included) in your ROBI COMB PRO's battery housing. Use of alkaline batteries is recommended. (See figure 1.)
2. This appliance suits any hair type.
3. Perform the treatment on clean, dry and oil free hair only.
4. Before the treatment, comb the hair thoroughly with a regular comb to free all knots.
5. Using the switch (see figure 2) operate the appliance in one of the two modes. The buzzing sound confirms that it is operating properly.
6. Hold the ROBI COMB PRO at a 45 degree angle and start combing with the combteeth near the hair-roots and the switch pointing upwards. (See figure 3.)

In mode 1: (identification) the appliance buzzes until it identifies and destroys a louse: at this point the buzz stops. Clean the comb with the cleaning brush. As soon as the combing unit is cleaned, the buzz recurs and you may continue. In this mode, you can quickly and easily identify the presence of lice.

In mode 2: (elimination) the appliance buzzes and destroys lice. After destroying the first louse, the buzzing tone changes while lice elimination continues. When the power supply runs out, the buzz stops. At this point clean the combing unit thoroughly, using the cleaning brush (figure 5). The buzz will return only when the unit is completely clean and free of lice, dandruff or any other dirt (figure 5).

7. In the event that the buzzing tone does not return, remove the comb housing by pressing the push buttons on each side of the appliance (figure 6.) Remove the combing unit from the housing (figure 7) and clean it thoroughly, using the combing brush. You may moisten the brush with alcohol (medical alcohol only). Once clean, dry the combing unit and brush thoroughly before re-use. A hair dryer can be used to reduce drying time.

8. The ROBI COMB PRO is highly effective in eliminating adult lice. Therefore, we recommend repeating treatments at least every day or two, for two weeks. You thus ensure elimination of small young lice, which escaped previous treatments, before they reach maturity and lay new eggs.

9. Every time your child complains of itching, we recommend combing his/her head with ROBI COMB PRO, close to the scalp. There is no need to run the comb over the entire hair-length.

10. If you will use ROBI COMB PRO for more than one person, then it should be cleaned thoroughly between each treatment. Follow the instructions under point 7.

11. Liceguard ROBI COMB PRO is especially efficient for grown lice and large nymphs. This means, that it is important to continue the treatment, ever other day until the remaining eggs have hatched. This will normally take 10-14 days. If ROBI COMB PRO is used according to the instructions, the lice will never become mature to reproduce. By using ROBI COMB PRO regularly, you will quickly detect the new lice. It is often the case, that children will have lice several times over a short period, since many parents are not consistent in combing their children in the necessary time frame.

ROBI COMB PRO makes it much easier to eliminate the problem.

figure 1

figure 2

figure 3

figure 4


figure 5

figure 6

figure 7

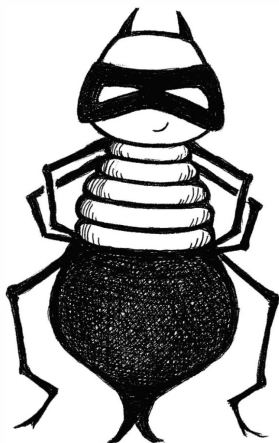FIND THE BANDITS!